

# FLAEM<sup>®</sup>

A HEALTHY LIFESTYLE

# AIRMIST

Aerosoltherapy unit for intermittent seasonal

## PRODUCT DESCRIPTION

- Piston compressor driver nebuliser
- Flaem's patented system guarantees comfortable usage
- Silent and rapid
- Effective treatment of the entire respiratory tract
- It can be used to administer any medication for aerosoltherapy
- Pneumatic circuit and accessories manufactured from biocompatible and anallergic materials

## EQUIPPED WITH

- Nebuliser RF6 basic with non-spill system, in vitro characterized by TÜV
- Connection hose
- Nebulisation manual control
- Mouthpiece
- Child Nasal prong
- Adult Nasal prong
- Adult mask in soft medical PVC



**!** Before switching the unit on, carefully read the "Instructions for Use" manual.

**F** Avant d'allumer l'appareil, consulter soigneusement le mode d'emploi.

**NL** De gebruiksaanwijzing aandachtig lezen voor het opstarten van het apparaat.

**D** Vor Einschalten des Inhalationsgerätes sorgfältig die Bedienungsanleitung lesen.

**E** Antes de encender el aparato, lea detenidamente el manual de instrucciones.

**RUS** Перед запуском машины, необходимо внимательно ознакомиться с руководством «Инструкции по пользованию».

**AR** استخدام الجهاز، على الكتيب "تعليمات الاستخدام".

**NEW**  
Nebuliser with non-spill  
system, in vitro  
characterized by  
**TÜV**



Convenient carrying handle



Convenient accessory compartment





Air purified by professional filter

**MADE IN ITALY**

[www.flaemnuova.it](http://www.flaemnuova.it)

## TECHNICAL SPECIFICATIONS


Power supply/Power:	<b>230 V ~ 50 Hz/140 VA</b> Other voltages are available upon request
Reciprocating piston compressor (designed for use without lubrication and maintenance)	
Max pressure:	<b>2.5 bar (approx.)</b>
Compressor air output:	<b>10 l/min (approx.)</b>
Sound level (at 1 m):	<b>58 dB (A) (approx.)</b>
Dimensions:	<b>27 (length) x 22 (width) x 11 (height) cm.</b>
Weight:	<b>2.200 Kg</b>
External air inlet and with a replaceable high efficiency filter (50/70 µm)	
In compliance with Dir. 93/42/EEC (for the 230 V ~ 50 Hz unit only):	 0051
Approvals:	
Warranty:	<b>Please refer to your local dealer</b>

Can be used with....




## JET NEBULISER TECHNICAL FEATURES

<b>RF6 BASIC nebuliser</b>	Average values measured with 0.9% saline solution in compliance with Flaem internal procedure. Certified with TSI Aerosizer 3220 computerised laser system in compliance with Flaem internal procedure.	In vitro characterization carried out at the Inamed Research GmbH & Co. KG on behalf of TÜV Rheinland Product Safety GmbH - Germany in compliance with the new European Standard for aerosol therapy appliances, Standard EN 13544-1. Further details are available upon request.
<b>Minimum capacity</b>	<b>2 ml</b>	
<b>Maximum capacity</b>	<b>8 ml</b>	
<b>Flow to nozzle</b>	<b>6.2 l/min</b>	
<b>Operating pressure</b>	<b>0.9 bar</b>	
<b>Max delivery</b>	<b>0.42 ml/min (approx)</b>	
<b>MMAD</b>	<b>1.53 µm</b>	<b>3.5 µm</b>
<b>Particles % &lt; 5 µm</b>	<b>99%</b>	<b>69%</b>

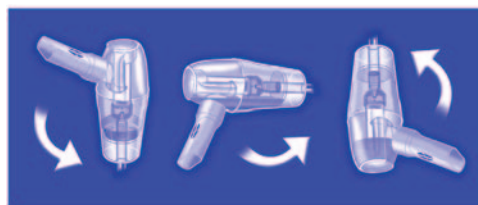

RF6 basic is also suitable for pediatric use, since the non-spill system prevents medication from spilling out even if the nebuliser is accidentally turned upside-down. (Efficacy tested up to 5 ml, the average dose recommended for normal use).

For further information on aerosol therapy and Flaem Nuova products please visit our website [www.flaemnuova.it](http://www.flaemnuova.it)

<b>AIRMIST F700 CLASSIC Series</b>	Code	<b>A009P00</b>	 8 018294 005238
	Master carton	<b>6 piece</b>	

Manufacturer:

**FLAEM NUOVA S.p.A**  
Elettromedicali/Elettrodomestici  
Via Colli Storici, 221-223-225  
25010 S. Martino della Battaglia  
Brescia (ITALY)  
Phone (39) 030 9910168  
Fax (39) 030 9910287  
E-mail: [info@flaemnuova.it](mailto:info@flaemnuova.it)  
[www.flaemnuova.it](http://www.flaemnuova.it)



Company with certified quality system

Export Department of FLAEM NUOVA



3, Via del Passero - 20147 Milano - Italy  
Phone (39) 02 4120769  
Fax (39) 02 48376295  
E-mail: [info@medexpointernational.com](mailto:info@medexpointernational.com)